



Premier Medical Necessity (PMN) (8.8)

Please read entire document before installing

System and software prerequisites

- You must be on v8.7 of PMN before installing this release.
- Your i5/OS must be at V6R1M0 or higher for this release to function properly.

Any program found in your HS library that is also found in the release library will need to be manually updated from the release library.

Download/Upload of save file and restore library

The save file for PMN v8.8 can be downloaded from our website under the Medical Necessity Product. www.Premsoft.com

The zipped file downloaded from the web page will need to be unzipped and the save file uploaded to your System i5 (iSeries) after which the release library will need to be restored. There are several methods you can use to accomplish this. Below is an example of using FTP to do so.

1. Extract the file pmn88.savf to a temporary directory. We'll use c:\temp in this example.
2. Create a save file on your System i5.
`CRTSAVF FILE(QGPL/PMN88)`
3. Open a DOS window and type "cd c:\temp"

4. In the same DOS window, use the 'ftp' command in Windows, to FTP 'pmn88.savf' into PMN88 in QGPL. Be sure to use the 'bin' command to change to binary mode. Here is a sample FTP session:

```
C:\temp>ftp nnn.nnn.nnn.nnn
Connected to nnn.nnn.nnn.nnn.
220-QTCP at (System i5 name).
220 Connection will close if idle more than 5 minutes.
User (nnn.nnn.nnn.nnn:(none)): userid
331 Enter password.
Password: XXXXXX
230 USERID logged on.
ftp> bin
200 Representation type is binary IMAGE.
ftp> put pmn88.savf QGPL/PMN88
200 PORT subcommand request successful.
150 Sending file to member PMN88 in file PMN88 in
library QGPL.
250 File transfer completed successfully.
8668704 bytes sent in 22.76seconds 380.19Kbytes/sec)
ftp> quit
221 QUIT subcommand received.
```

5. From a System i5 command line, run the following command.


```
RSTLIB SAVLIB(PMN_V8.8) DEV(*SAVF) SAVF(QGPL/PMN88)
```

6. On a command line, run the following command.

```
DLTF FILE(QGPL/PMN88)
```

Install Tasks

1. Sign on to your System i5 with a user profile that has at least SECADM authority.
2. Set your library list for the PMN environment. **Should include PH/MR/EC/SC if installed.**
3. Add the release library to your library list (PMN_V8.8)
4. On a command line type **\$_install** and press enter. The installation notice screen displays.

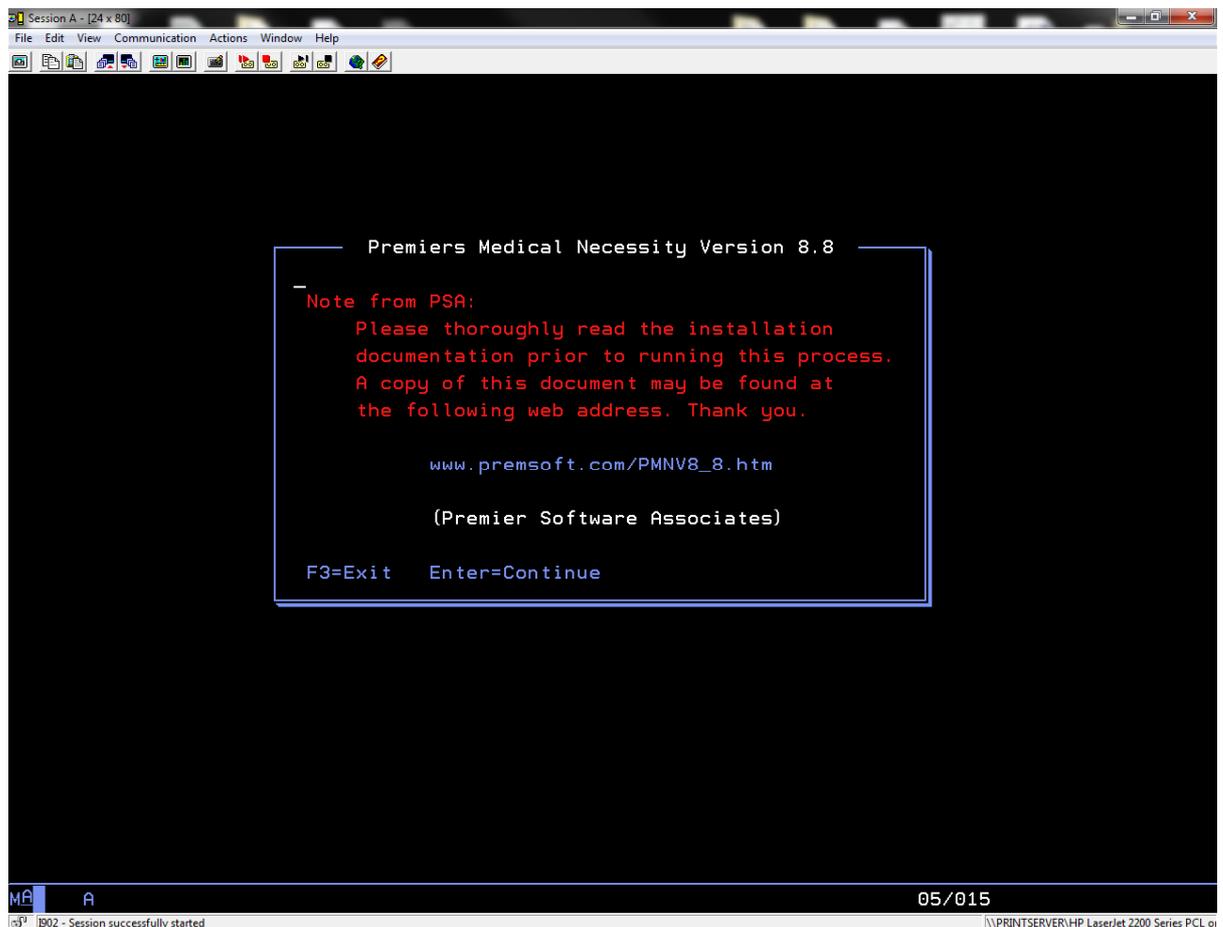


```
Session A - [24 x 80]
File Edit View Communication Actions Window Help
Premiers Medical Necessity Version 8.8
Note from PSA:
Please thoroughly read the installation
documentation prior to running this process.
A copy of this document may be found at
the following web address. Thank you.

www.premsoft.com/PMNV8_8.htm

(Premier Software Associates)

F3=Exit Enter=Continue
05/015
IB02 - Session successfully started \\PRINTSERVER\HP LaserJet 2200 Series PCL or
```


5. Press Enter to continue. The environment library screen displays

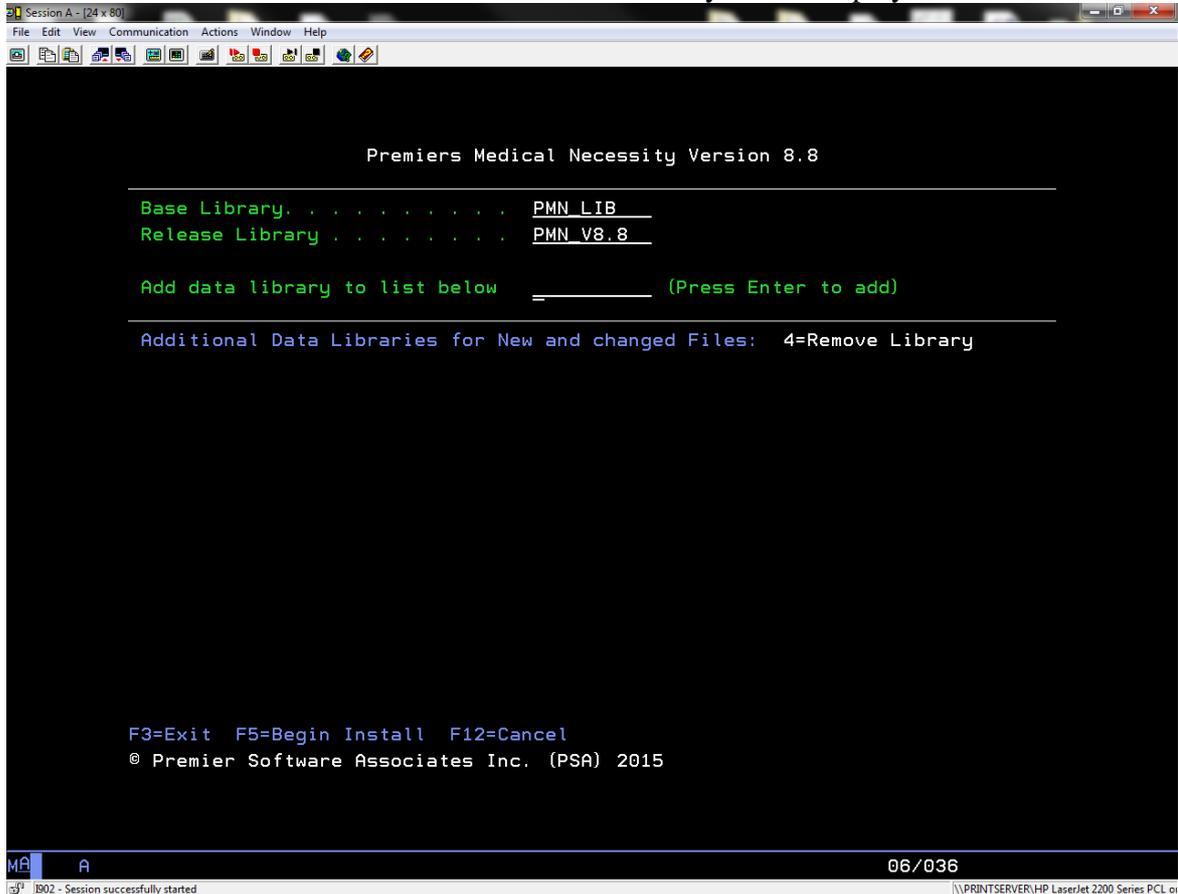


6. Fill out the fields on the screen according to the following table

Field	Description
Base Library	Enter the name of the library containing the current version of PMN which is to be updated to v8.8. This library contains the PMN programs and may contain data files in a single entity environment
Release Library	Enter the name of the release library, PMN_V8.8 (or equivalent if you restored the release under a different name).
Add data library...	For this particular install of Version 8.8 The files being changed do not go into additional data libraries. <u>DO NOT add any additional data libraries for this install. Thank you.</u>
	*Note – If additional libraries are added, then you will be presented with an additional screen listing them during the install.

7. Press **F5** to start the install process. The installation verification window will display.

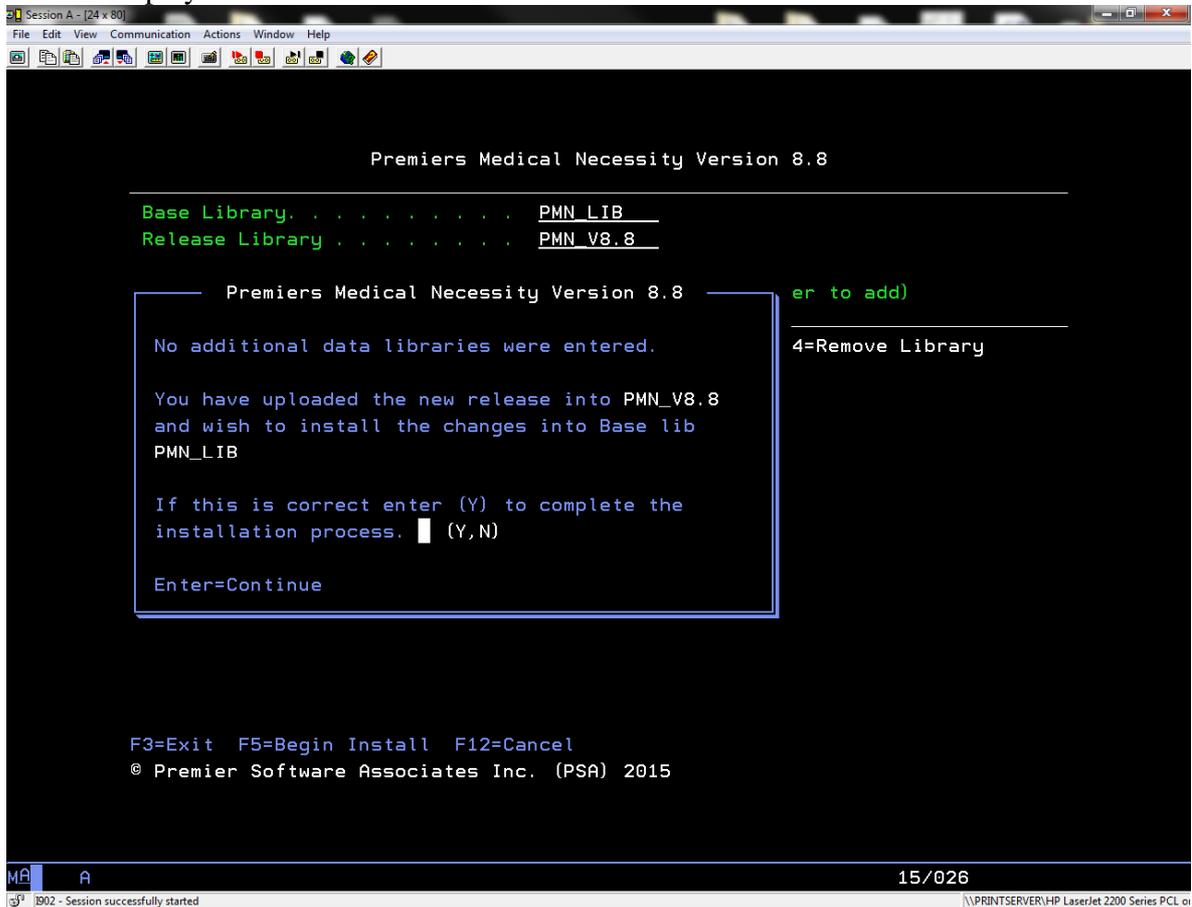

8. Enter a "Y" and press enter to install the release.
9. When the installation is complete the following window displays.

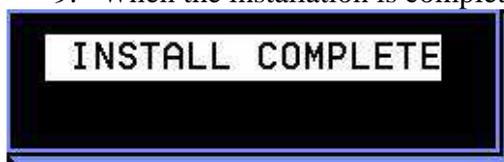

10. Press enter to exit the installation program.

11. Query to review the file \$_INSTERP in PMN_V8.8. It will contain records of any errors which might have occurred during the install.
 - a. The file will list any commands that the install program ran which were not able to complete successfully.
 - b. Depending on your environment, you may receive the following errors on the listing. These are normal and do not require any remediation.

CRTRPGPGM (IN0003R)

- c. Any other commands listed in this file may require remediation. Analyze your environment to determine the reason for the failing command.
12. Repeat steps 1-11 for any additional PMN environments which do not share the base library just upgraded.
13. **Additional Notes.**

Additional Notes

The following programs will not compile if the current system does not have the MS4 scheduling software installed.

- a. IN0003R
- b. This install will clear the ICD10 file in PMN. You will need to reload this file by going to the Medical Necessity Main Menu and following these steps.
 1. Take option 1 (Maintenance Menu)
 2. Take option 22 (File Maintenance Menu #2)
 3. Take option 6 (Synchronize ICD10 File)
 4. Take Option 5 (Rebuild ICD10 Search Index)
 5. Now option 4 (ICD10 Maintenance) will be loaded.

Any V8.8 programs existing in HS libraries will need to be recompiled.